
GlucoMemory Software Schritt für Schritt Step by step Passo dopo passo

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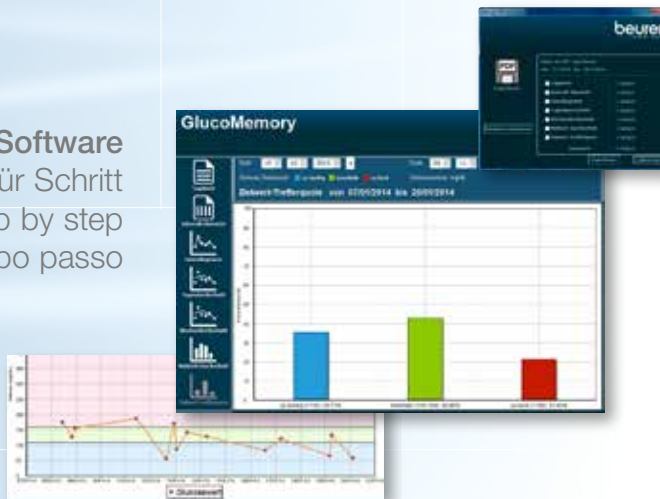Contents: GlucoMemory software

1 Overview	30
1.1 This manual	30
1.2 Notes on reading the manual	30
1.3 Intended use of GlucoMemory software	30
2 Introduction to GlucoMemory	31
2.1 GlucoMemory software	31
2.2 Compatible blood glucose monitors	31
2.3 System requirements	31
2.4 User interface operation	32
2.4.1 Components of the user interface	32
2.4.2 Buttons for selecting the evaluation	32
2.4.3 Administration functions	33
2.4.4 Date selection	34
2.5 Colours and backgrounds	35
3 Overview of software versions	36
3.1 Embedded version (GL50, GL50 evo only)	36
3.2 PC version	36
4 Installation and initial use of the software	36
4.1 Installing the PC version	36
4.1.1 Installation directly on the monitor (GL50, GL50 evo only)	37
4.1.2 Installation using the PC kit (GL44, GL40)	37
4.1.3 Installation using the PC kit (GL34, GL32)	37
4.2 Putting the embedded version into operation (GL50, GL50 evo only)	37
4.3 Connecting the monitor to the PC and starting the software	38
4.4 Downloading the blood glucose measurements from the monitor (PC version)	38
4.5 Disconnect the monitor from the PC	38
4.6 Uninstalling the PC version	39
5 Setting diabetes management	39
5.1 Setting blood glucose units	40
5.2 Setting blood glucose targets	41
5.3 Setting time blocks for before and after meals	42
6 Evaluating blood glucose data	43
6.1 List of blood glucose data: log book	43
6.1.1 Displaying data in the log book	43
6.1.2 Editing data in the log book	44
6.1.3 Removing data from the log book	45
6.2 Evaluation: range overview	46
6.3 Evaluation: line chart	47
6.4 Evaluation: average day trend	48

6.5 Evaluation: average week trend.....	49
6.6 Evaluation: avg meal histogram.....	50
6.7 Evaluation: target ratio.....	51
7.0 Managing blood glucose data.....	52
7.1 Exporting data.....	52
7.1.1 Exporting to PDF file.....	53
7.1.2 Exporting to CSV file.....	53
7.2 Print out data on the printer.....	53
8 Error messages and their solution.....	54

1 OVERVIEW

1.1 This manual

This manual provides a thorough description of the operation and configuration of **GlucoMemory** software. It includes information on data exchange and export as well as the creation of evaluations. The manual is intended for all users, GPs and healthcare professionals who work with **GlucoMemory** software.

1.2 Notes on reading the manual

To make reading the manual easier, the following typographic aids are used:

Aid	Description
Proper name	Product names, the designation of buttons and symbols, texts in dialogues and emphasised terms appear in bold.
Input text	Texts that have to be input by the user are displayed in this special character format.
<i>Italic</i>	File names are written in italic font.
 Note	Notes and tips that make work easier.
 Caution!	Safety note on serious potential errors.

1.3 Intended use of GlucoMemory software

The **GlucoMemory** software is primarily for the manual, computer-supported documentation of metabolic data. The principal purpose is as a contemporary replacement for the blood pressure diaries that were previously handwritten by patients. The software is intended exclusively for persons who are well informed about type 2 diabetes and the associated treatment strategies. In case of doubt, assistance should be sought from a doctor who can advise the patient about how to use the software to treat diabetes.

GlucoMemory should only be used in the context of the above intended use. In case of any uncertainty about the diabetes treatment program, advice must be sought from a doctor. Therapeutic decisions and changes should only be made in consultation with a doctor. The GlucoMemory software is intended exclusively for the supplementary documentation, presentation and archiving of therapeutic data.

Caution!

Statistics, analyses and evaluations produced by GlucoMemory should never be used as the sole basis for therapeutic decisions, but should be viewed only as a supplementary, digital blood pressure diary. The use of GlucoMemory releases neither the user nor the doctor from their responsibility to verify that the data obtained from the program is correct.

2 INTRODUCTION TO GLUCOMEMORY

2.1 GlucoMemory software

GlucoMemory is an advanced software for evaluating your blood glucose measurements. The evaluation software enables you and your GP or healthcare professional to better monitor changes in your blood glucose level and manage your diabetes.

The **GlucoMemory** software enables you to view your blood glucose values on the PC and offers a variety of evaluation options. A blood glucose reading contains the date, time and the blood glucose value recorded. **GlucoMemory** informs you of your blood glucoserelated values, displays these in lists, tables and charts, thus supporting your diabetes management.

2.2 Compatible blood glucose monitors

Beurer GL50 evo, GL50, GL44, GL43, GL40, GL34 and GL32 blood glucose monitors.

2.3 System requirements

In order to work with the **GlucoMemory** software, you will need a standard PC that fulfils the following minimum requirements:

Operating system	<ul style="list-style-type: none">• Windows: Microsoft® Windows® 7, Vista (SP1), XP (SP3)• Macintosh: Windows® emulation on an Apple Mac computer, e.g. Fusion, Parallels or Virtual PC
PC	<ul style="list-style-type: none">• Prozessor: Intel Pentium II 300 MHz CPU or higher• Main memory: 32 megabytes (MB) or greater
Interfaces	1 USB interface

The following resolutions are recommended for optimal use of **GlucoMemory**:


- Screen resolution: 1024 x 768 pixels
- The best resolution is 96 DPI (dots per inch).

2.4 User interface operation

This chapter describes the operation of individual components and operating elements of the GlucoMemory user interface.

2.4.1 Components of the user interface

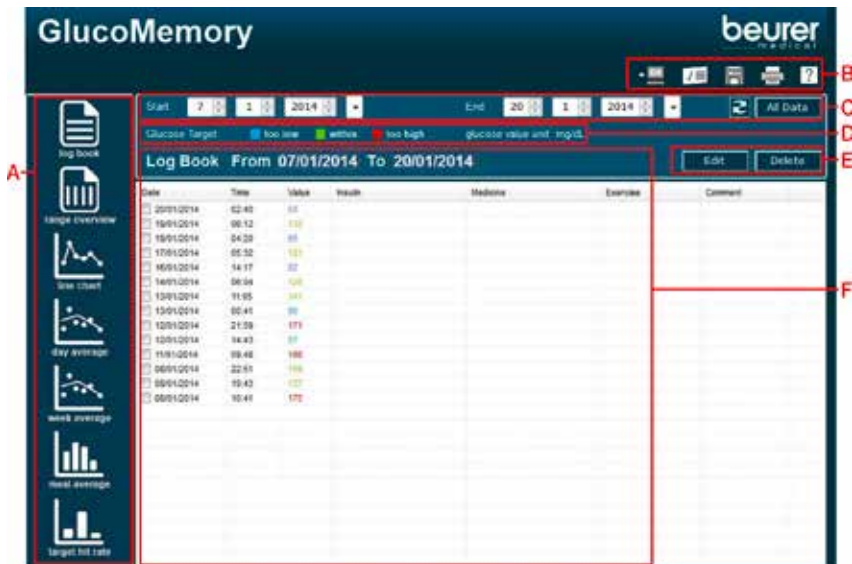The work area is made up of the following components:

- A** Buttons for selecting an evaluation
- B** Toolbar with setting and administration functions as well as the Download button (only in PC version)
- C** Date selection
- D** Display of target colours and measuring unit (mg/dL or mmol/L)
- E** Buttons for editing and deleting a measurement
- F** Display area for the selected evaluation

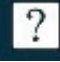
2.4.2 Buttons for selecting the evaluation

Open an evaluation by clicking the corresponding button. The evaluation is then shown in the display area. The button for the selected evaluation appears in grey.

 log book	log book	List of the blood glucose values sorted by date and time (the most recent value appears at the top of the list)
 range overview	range overview	List of the blood glucose values sorted by time blocks before and after meals
 line chart	line chart	Line chart showing the daily blood glucose values
 day average	average day trend	Chart showing the daily blood glucose values (plotted as dots) and the average values for each time block (plotted as a line)
 week average	average week trend	Chart showing the daily blood glucose values (plotted as dots) and the average values for each week (plotted as a line)
 meal average	avg meal histogram	Bar chart showing the average values for blood glucose before and after each meal
 target hit rate	target ratio	Bar chart showing the percentage shares according to the blood glucose target (too low, within, too high)

2.4.3 Administration functions

Call up an administration function from the toolbar by clicking it.

	Analysis Setting	Set blood glucose targets, time blocks and glucose unit for blood glucose management
	Export	Export data to a PDF or CSV file
	Print	Print out data on a printer
	Help	Manufacturer information, software version, support information and manual

	Download	Only in PC version: Transfer the blood glucose readings from the blood glucose monitor to the PC
---	----------	---

2.4.4 Date selection

Using date selection, you can define a time block for the values of an evaluation. Date selection is supported by all evaluation types and can be set differently for each evaluation type. The components of date selection are:

- A** Start date (with counter fields for day/month/year)
- B** End date (with counter fields for day/month/year)
- C** Refresh the evaluation in accordance with date range
- D** Display all stored data in the evaluation

Note

The [drop-down]  symbol beside the start and end date opens a calendar for convenient date selection.

Selecting date range and refreshing the evaluation (for example: January, 2014)

1	Enter the start date in the Start fields or select this date using the calendar (for example: 14 1 2014).
2	Enter the end date in the End fields or select this date using the calendar (for example: 22 1 2014).
3	Click Refresh  . The evaluation that was just opened is displayed with the values for the updated date range. The new time block is displayed on the right next to the title of the evaluation.

Displaying all data in an evaluation

1

Click **All Data** . The evaluation that was just opened takes into account all stored data. The new time block is displayed on the right next to the title of the evaluation.

2.5 Colours and backgrounds

The evaluations use various colours and backgrounds to present the blood glucose data.

Evaluations: log book and range overview

In these evaluations, the software uses three different colours to classify the data.

Target colour	Blood glucose measurement	Example
	> 160 Above the blood glucose target	
	110 - 160 Within the blood glucose target	
	< 110 Below the blood glucose target	

 **Note:** The target can be set individually using **Settings** .

Evaluations: line chart, average day trend and average week trend

In these evaluations, the software uses different coloured backgrounds to illustrate the data range.

Background	Blood glucose measurement	Example
	> 160 Above the blood glucose target	
	110 - 160 Within the blood glucose target	
	< 110 Below the blood glucose target	

 **Note:** The target can be set individually using **Settings** .

3 OVERVIEW OF SOFTWARE VERSIONS

The software is available in different versions: Embedded version and PC version.

ⓘ Note: The language of the GlucoMemory user interface is determined by the language of the operating system. This applies for the embedded version as well as the PC version. The languages available are English and German.

3.1 Embedded version (GL50, GL50 evo only)

This refers to the software integrated in the blood glucose monitor (embedded version). The software runs directly on the monitor. To avoid the spreading of PC viruses, no data can be saved on the monitor. This means that if you make changes to the embedded version, such as adding or deleting data, these changes are not saved on the monitor when the software is closed.

ⓘ Note: Once the measurement memory on your blood glucose monitor is full, new measurements automatically replace the oldest values stored. By installing the PC version (see chapter „3.2 PC version“), you can save all your values for the long term on the PC and add further details.

3.2 PC version

In addition to the embedded version, there is a version that can be installed on the PC.

- With the GL50, GL50 evo blood glucose monitors, this PC version is also directly on the monitor.
- With the GL44, GL40, GL34 and GL32 blood glucose monitors, the PC version is included as a CD ROM in the PC kit (optionally available). You can download the software also free of charge under www.beurer.com.

You can transfer your measurements to the PC version of the software. In the PC version, you can permanently save changes, such as adding or deleting data or adjusting the settings.

ⓘ Note: Please note that values deleted from the software are not deleted from the monitor!

4 INSTALLATION AND INITIAL USE OF THE SOFTWARE

Caution!

- Read this chapter carefully before connecting the monitor to the PC. Follow the instructions relating to installation and initial use of the software very closely and in the correct order.
- Check the battery status of the monitor. The monitor must contain a functional battery.

4.1 Installing the PC version

ⓘ Notes:

- When the blood glucose values are transferred from the monitor to the software, the existing data is not overwritten. The new data is added to the software.
- When you close the **GlucoMemory** software, all the data is backed up in the software by confirming the security prompt.

4.1.1 Installation directly on the monitor (GL50, GL50 evo only)

1	Install the PC version in the file explorer using the executable file <i>GlucoseMemory_Setup.exe</i> .
2	Follow the steps in the GlucoseMemory Setup Assistant .

4.1.2 Installation using the PC kit (GL44, GL40)

1	Insert the supplied CD ROM in the CD drive. Install the PC version in the file explorer using the executable file <i>GlucoseMemory_Setup.exe</i> . The GlucoseMemory Setup Assistant opens automatically if the Auto-Run function is supported by your computer.
2	Follow the steps in the GlucoseMemory Setup Assistant .

4.1.3 Installation using the PC kit (GL34, GL32)

1	Insert the supplied CD ROM in the CD drive. Install the PC version in the file explorer using the executable file <i>GlucoseMemory_Setup.exe</i> . The GlucoseMemory Setup Assistant opens automatically if the Auto-Run function is supported by your computer.
2	Follow the steps in the GlucoseMemory Setup Assistant .


That your blood glucose monitoring devices GL34 und GL32 can communicate with your PC, you need to install a driver file.

1	Start the file "CP210x_VCP_Win_XP_S2K3_Vista_7.exe" in the folder „driver GL34_GL32_BGL60“.
2	Follow the steps in the Silicon Laboratories Driver Setup Assistant .

4.2 Putting the embedded version into operation (GL50, GL50 evo only)

The software is integrated in the Beurer GL50, GL50 evo blood glucose monitors and requires no additional installation of system components on local PCs. For further steps, see next chapter "4.3 Connecting the monitor to the PC and starting the software" 4.3.

4.3 Connecting the monitor to the PC and starting the software

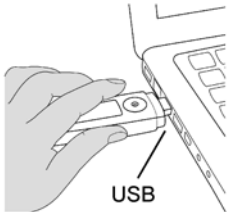1	Ensure that no other blood glucose monitor is connected to the PC. Only one monitor can be connected to the PC at a time.
2	<p>Connect the blood glucose monitor to the PC via the USB interface as described in the monitor's Instructions for use.</p> <p>The image shows the USB connection with the GL50 monitor as an example.</p> 
3	<p>Depending on the software version, proceed as follows:</p> <p>PC version: Start the software via the start menu or open the software using the desktop icon.</p> <p>Embedded version (GL50, GL50 evo): Depending on whether or not the computer supports the Auto-Run function, two situations can occur:</p> <ul style="list-style-type: none">• If Auto-Run is activated, the integrated GlucoMemory software automatically starts and displays the log book with all data as the start screen.• If Auto-Run is deactivated, you must start the GlucoMemory manually. To do so, open the Windows file explorer, click <i>GlucoMemory</i> and then the executable file <i>GlucoMemory.exe</i> in the subfolder \ <i>GlucoMemory</i>. The software displays the log book with all data as the start screen.

4.4 Downloading the blood glucose measurements from the monitor (PC version)

1	Connect your monitor and start the software as described above.
2	To transfer your new measurements to the software, click the Download  symbol in the toolbar.

4.5 Disconnect the monitor from the PC


To disconnect the monitor from the PC properly, proceed as follows.

1	<p>Close the GlucoMemory software using the Close  symbol (in the top, right corner of the window). Confirm the security prompt.</p> <p>Note on embedded version (GL50, GL50 evo)</p> <p>If you have added or deleted measurements in the GlucoMemory software, it is advisable to export these before closing (see chapter "7 Managing blood glucose data"). The changes are not saved when the software is closed.</p>
---	--

2	Remove the monitor using the Windows function “Safely Remove Hardware”  in the taskbar along the bottom edge of the screen.
3	Remove the monitor from the USB interface of the PC or docking station. <i>GlucoMemory</i> then no longer appears under <i>Desktop</i> in the file explorer.

4.6 Uninstalling the PC version.

You can remove the PC version from your PC at any time.

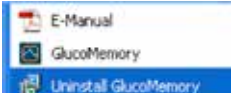1	Call-up: ⇒ Start menu ⇒ Programs ⇒ GlucoMemory ⇒ Uninstall GlucoMemory 
2	Confirm the security prompt with “Yes”. 

5 SETTING DIABETES MANAGEMENT

Notes:

- With the embedded version (GL50, GL50 evo only), it is not possible to save. Your changes are not saved on the monitor. Perform a data export (see chapter “7 Managing blood glucose data”) and/or install the PC version to retain your changes.
- The language of the GlucoMemory user interface is determined by the language of the operating system. This applies for the embedded version as well as the PC version. The languages available are English and German.

5.1 Setting blood glucose units

Caution !

You must set the unit for the evaluation to match that of your blood glucose monitor.

1	Click the Settings  symbol in the toolbar. The Settings dialogue opens.
2	Change the selection in the Unit block if necessary. 
3	Save by clicking Save.

5.2 Setting blood glucose targets

You have the option to set the blood glucose targets. The targets are used in the various evaluations.

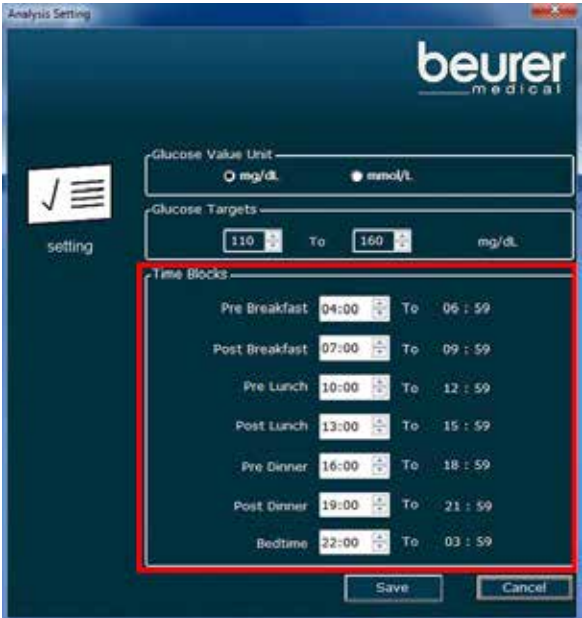
⚠ Caution!

- If you are unsure about the targets, ask your GP or healthcare professional.
- When the values are displayed on the PC, they are already linked to the preset targets. If you change the targets, the colour representations in the **GlucoMemory** software are adjusted.

1	Click the Settings  symbol in the toolbar. The Settings dialogue opens.
2	<p>The software provides standard values. If you wish, you can change the lower and/or upper value of your blood glucose target in the Glucose Targets block.</p> <p>⚠ Caution ! Be sure to note whether the unit mg/dL or mmol/L is selected!</p> 
3	Save by clicking Save.

5.3 Setting time blocks for before and after meals

You can define 7 time blocks that are before and after the various mealtimes of a day. The time blocks affect the evaluations for range overview (see chapter “6.2 Evaluation: range overview”), average day trend (see chapter “6.4 Evaluation: average day trend”) and avg meal histogram (see chapter “6.6 Evaluation: avg meal histogram”).

1	Click the Settings  symbol in the toolbar. The Settings dialogue opens.
2	<p>The software provides standard values for the time blocks. If you wish, you can change the times for the various meals according to your daily routine in the Time Block area. Each of the times indicate the start of the time block.</p> 
3	Save by clicking Save.

6 EVALUATING BLOOD GLUCOSE DATA


6.1 List of blood glucose data: log book

6.1.1 Displaying data in the log book

This list shows the detailed blood glucose readings, sorted by date and time. You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”). You can edit and delete the data as described in the following sections.

Call-up:

⇒ Automatically when opening the **GlucoMemory** software. Or click the **log book** button.

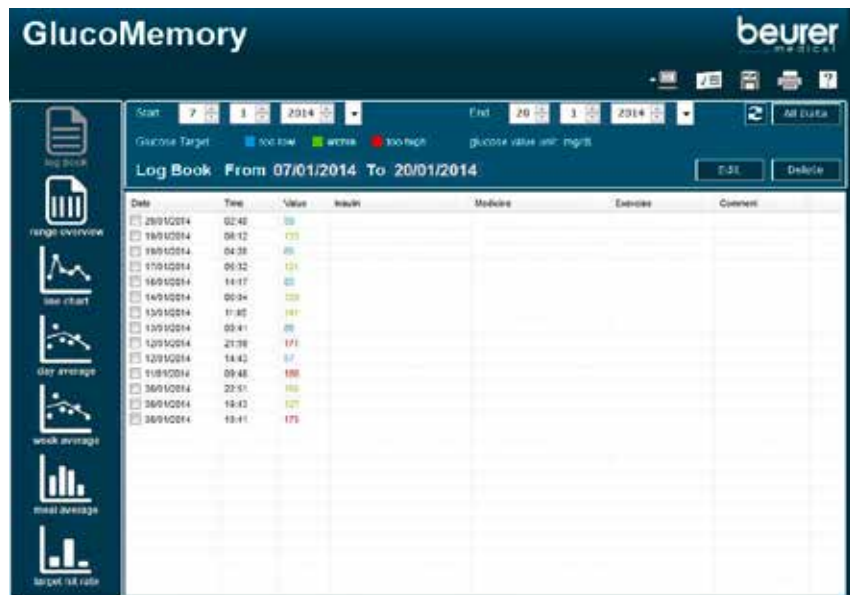
6.1.2 Editing data in the log book

You can add to blood glucose readings from the list. The insulin types and their dose, the medicine taken, the physical exercise and the comments can be changed. You can add one or several values at the same time.

- 1 You can edit one or several values at the same time:
 - To edit individual values: double-click the value.
 - To edit several values at the same time: mark the desired value in the list by clicking the boxes to the left of the blood glucose value. Then click Edit .

- 2 The **Edit** dialogue opens. Make the desired entries.

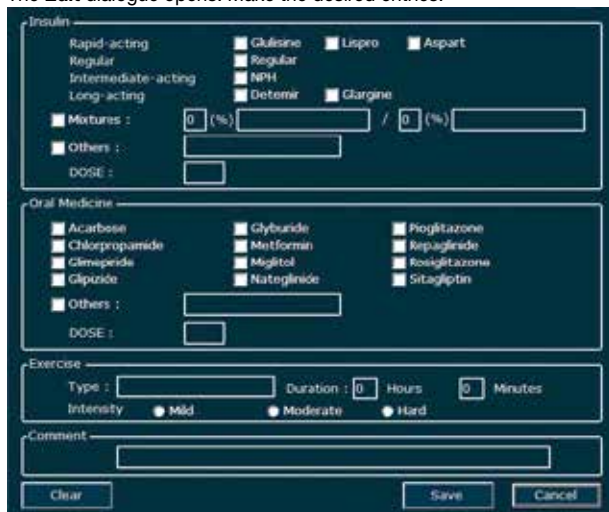
- 3 You can cancel your entries at any time using **Clear**. Then save your information by clicking Save. The entries are displayed in the relevant columns of the list.

<input type="checkbox"/>	14/01/2014	06:04	120		
<input type="checkbox"/>	13/01/2014	11:05	141		Metformin
<input type="checkbox"/>	13/01/2014	00:41	86		Metformin
<input type="checkbox"/>	12/01/2014	21:59	171		
<input type="checkbox"/>	12/01/2014	14:43	57		
<input type="checkbox"/>	11/01/2014	09:48	186		

6.1.3 Removing data from the log book

You have the option to delete one or several values at the same time. The deleted values are temporarily removed from the GlucoMemory software. However, the data continues to be saved in the monitor. If you would like to permanently delete the values, you must also delete them directly on the monitor using the Delete function (see description on deleting all data or individual measurements in the instructions for your blood glucose monitor).

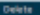1	Mark one or several values in the list by clicking the boxes to the left of the blood glucose value.
2	Click Delete  . Confirm the security prompt with OK. The list is then refreshed and the deleted values are no longer displayed.
3	If you wish to permanently delete the values, you must also remove these from your monitor.  Caution ! Always disconnect the monitor from the PC properly (see chapter “4.5 Disconnect the monitor from the PC”).


6.2 Evaluation: range overview

This evaluation shows the measurements with their allocation to the 7 time blocks before and after meals. You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

For information on how to set the times of the 7 time blocks individually, see chapter “5.3 Setting time blocks for before and after meals”.

Call-up:

⇒ Click the **range overview** button.

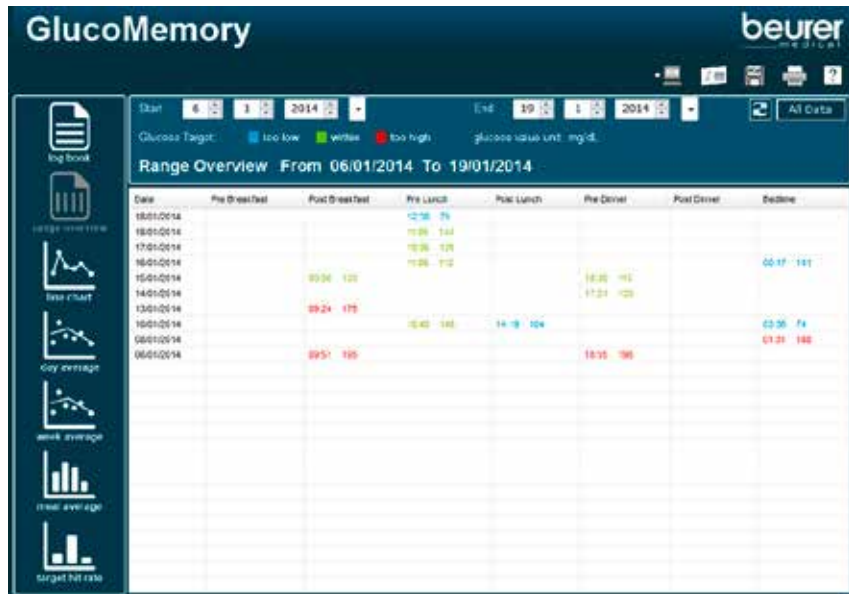6.3 Evaluation: line chart

This evaluation shows the course of the daily measurements in a line chart. You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

The blood glucose values are displayed in the chart as red dots. Each value is connected to the next with an orange line. The target areas are also filled with the coloured backgrounds (see chapter “2.5 Colours and backgrounds”).

Call-up:

⇒ Click the **line chart** button.

6.4 Evaluation: average day trend

This evaluation shows the average blood glucose value for all 7 time blocks (from left to right in the chart). You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

The blood glucose values are displayed in the chart as red dots. The average value is shown by an orange line. If a red dot, i.e. an individual blood glucose value, lies above or below the line, this indicates a high or low value compared with the average. The different target areas are also filled with the coloured backgrounds (see chapter “2.5 Colours and backgrounds”).

Call-up:

⇒ Click the **day average** button.

6.5 Evaluation: average week trend

This evaluation shows the average blood glucose value for the days of a week (from left to right in the chart). You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

The blood glucose values are displayed in the chart as red dots. The average value is shown by an orange line. If a red dot (individual blood glucose value) lies above or below the line, this indicates a high or low value compared with the average. The target areas are also filled with the coloured backgrounds (see chapter “2.5 Colours and backgrounds”).

Call-up:

⇒ Click the **week average** button.

6.6 Evaluation: avg meal histogram

This evaluation shows the average blood glucose value for various mealtimes (from left to right in the chart: breakfast, lunch, dinner and all mealtimes). For each mealtime, the average values before and after the meals are displayed side-by-side so that a direct comparison can be made.


The colour of the bars indicate the following:

- **Magenta:** average values before meals
- **Orange:** average values after meals

You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

Call-up:

⇒ Click the **meal average** button.

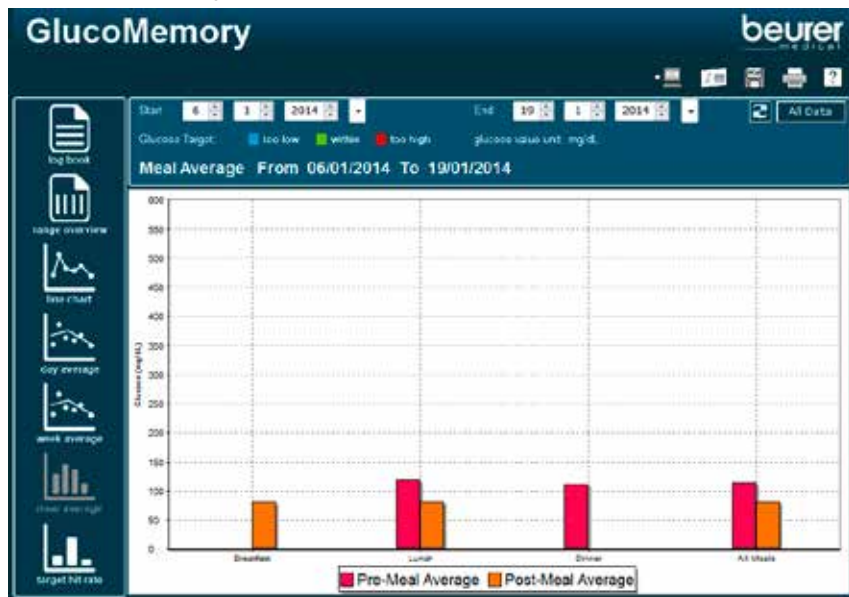6.7 Evaluation: target ratio


This evaluation shows the distribution of the blood glucose values compared with the defined target (see chapter “5.2 Setting blood glucose targets”).

- **Blue:** values below the target
- **Green:** values within the target
- **Red:** values above the target

You can select data from a specific date range or display all of the data (see chapter “2.4.4 Date selection”).

Call-up:

⇒ Click the **target hit rate** button.

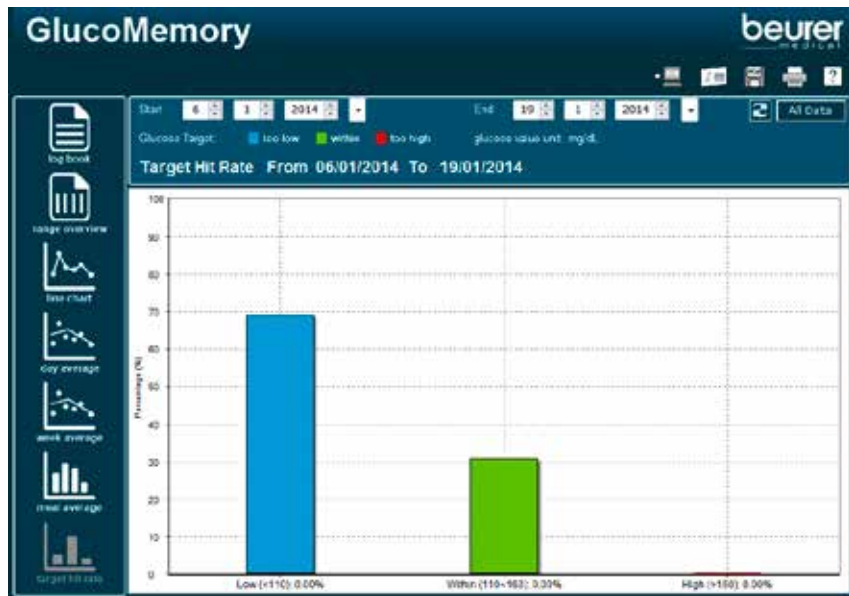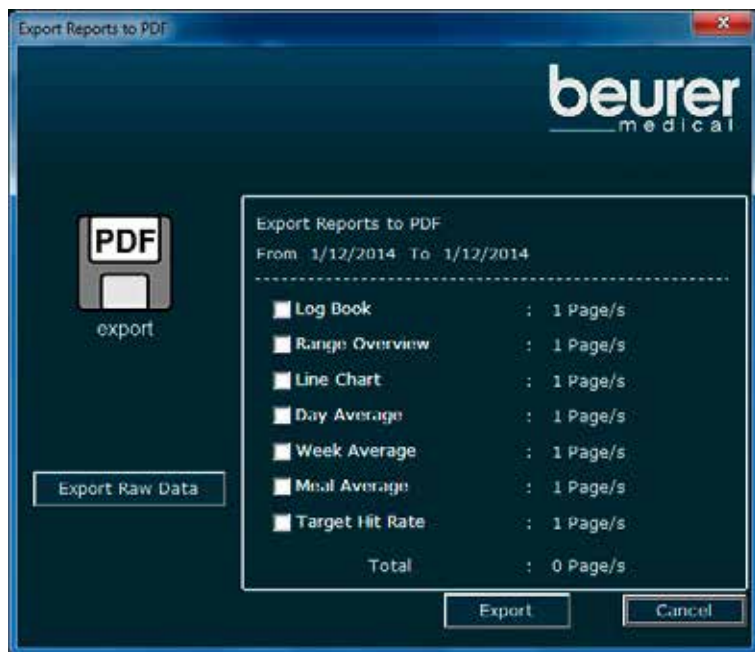
7.0 MANAGING BLOOD GLUCOSE DATA

To permanently backup and archive the data, you can export it to a PDF or CSV file (Microsoft Excel®) or print it out on a printer.

7.1 Exporting data

This chapter describes how you can export the data to a PDF or CSV file (Excel®).

1	In the date fields, select the desired time block for the data to be exported.
2	Click the Export  symbol in the toolbar. The Export Reports to PDF dialogue opens.
3	The evaluations are displayed in the PDF with the relevant number of pages. The total number of pages is shown below the individual pages.
4	Export to a PDF or CSV file: <ul style="list-style-type: none">• If you wish to export to a PDF file, please see chapter “7.1.1 Exporting to PDF file”.• If you wish to export to a CSV file (Excel®), please see chapter “7.1.2 Exporting”.

7.1.1 Exporting to PDF file

You can export all evaluations to a PDF file for the purposes of archiving and back them up in any memory area of a computer.

1	The Export Reports to PDF dialogue is opened. Activate the desired evaluations by clicking the boxes.
2	Click Export .
3	The file dialogue Save as opens. Enter the target folder and file names of the PDF file and click Save . The PDF file then appears in the target folder.

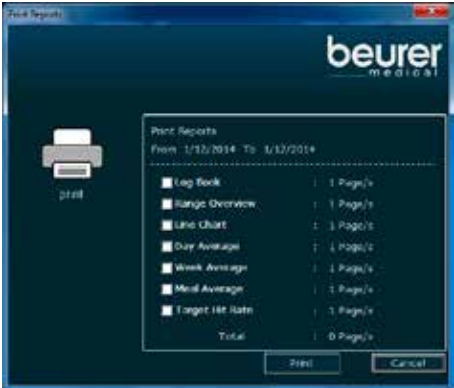
7.1.2 Exporting to CSV file

You have the option of exporting the data of a specific time block to a CSV file. This export firstly can be used for archiving your data and secondly enables you to perform further calculations and statistical evaluations in programs such as Microsoft Excel®.

1	The Export Reports to PDF dialogue is opened. Activate the desired evaluations by clicking the boxes.
2	Click Export Raw Data .
3	The file dialogue Save as opens. Enter the target folder and file names of the CSV file and click Save . The CSV file then appears in the target folder.

7.2 Print out data on the printer

You can print out all evaluations for the purposes of long-term archiving.

1	In the date fields, select the desired time block for the data to be printed.
2	Click the Print  symbol in the toolbar. The Print Reports dialogue opens. 

3	Activate the desired evaluations by clicking the boxes. Click Print .
4	In the Print dialogue, select the desired printer and if necessary the desired number of copies, then click OK . The evaluations are then printed out on the printer.

8 ERROR MESSAGES AND THEIR SOLUTION

Error message	Solution
USB connection failed.	The USB connection has failed. Reconnect the monitor to the PC (see chapter "4.3 Connecting the monitor to the PC and starting the software").
Edit item cannot contain any of the following characters: # * \$!	Ensure that you do not use the characters # * \$! when entering the data.
Time blocks exceed 24 hours.	This message appears if a time block exceeds 24 hours when setting the time blocks for mealtimes (see chapter "5.3 Setting time blocks for before and after meals"). Set the times according to your mealtimes.
Errors happened on period date setting.	During date selection, ensure that the end date is not BEFORE the start date. The end date must be later than the start date. For a description of date selection, see chapter "2.4.4 Date selection".

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